
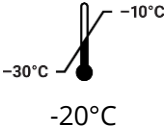








INSTRUCTIONS FOR SAMPLE COLLECTION AND SHIPMENT OF iCS-DIGITAL™ AND STR TESTS

Sample	Method	Volume & concentration	Temperature	
			Storage (optional)	Shipment
 DNA	<p>Harvest genomic DNA from the cell line to be analysed using an appropriate genomic extraction kit.</p> <p>We recommend the use of a column extraction kit for optimal digital PCR results (do not use quick extraction protocols).</p> <p>Determine the concentration and quality of genomic DNA samples using an appropriate method.</p>	<p>DNA concentration: ≥ 10 ng/μL if dosed by Qubit or ≥ 50 ng/μL if dosed by Nanodrop</p> <p>Volume: ≥ 40 μL*</p>		 Ambient temperature
 Cells	<p>Dissociate cells (with enzymes, EDTA, etc.), centrifuge them (1200 rpm; 5 min), and discard the supernatant. Wash the cell pellet with PBS, centrifuge and discard supernatant.</p> <p>Frozen cell pellets must be sent on dry ice.</p>	<p>Dry pellet ≥ 500,000 cells</p>		 Dry ice

*** Warning:** if the volume sent is **< 40μL**, it will be **levelled to 40μL with H2O** then **quantified with the Qubit fluorometer**.

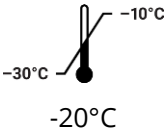
Depending on how many tests you run on the same cell line, please add 500,000 cells per additional test.

Fill in the [Sample release form](#) carefully for all the samples to be tested and [send it by e-mail in an Excel format to services@stemgenomics.com](#).

Shipment (The following instructions must be strictly followed)

Your sample(s) must be sent in four layers of packaging as follows:

- 1) Sample(s) must be in a [leakproof tube](#)
- 2) Place the tube(s) in a [first sealed bag with absorbent material](#) (for example paper towel). The maximum capacity allowed is 1L.
- 3) Place the first sealed bag in a [second bigger and sealed bag](#)

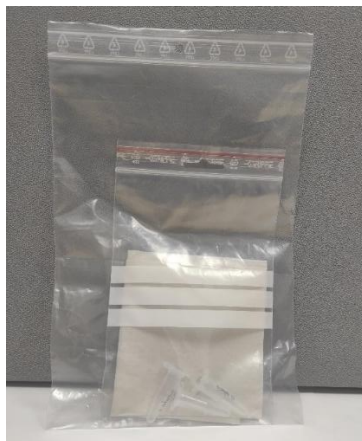



- 4) Place the sample in a rigid packaging or a bubble wrap envelop.

Send it to:

Stem Genomics - Testing Service
Cap Sigma
1682 rue de la Valsière
34790 Grabels
France
VAT n°: FR79841152077
EORI n°: FR841152077

Phone number and email address to share
with the transporter:

Phone: +33 420 900 201

Email: services@stemgenomics.com

Total declared value: 10 €

Description of Goods: samples, non-infectious, non-hazardous, not for medical purposes, for research purposes only, not for resale. Goods packed in compliance with all applicable regulations. Not for human use.

HS Code: 3002.59 (or 3002.5900)

The proforma invoice / commercial invoice must clearly indicate the mention "no commercial value" and "value for customs purposes only"

Cell pellets must be sent on dry ice and therefore must comply with Dry Ice Shipping Regulations. Please contact your carrier.